



Presto Stantest Pvt. Ltd.



I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India
P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com

- Faridabad • Sonipat • Kolkata • Mumbai • Pune • Ahmedabad
- Chennai • Bangalore • Hyderabad

www.prestogroup.com



Flammability Tester (UL-94)

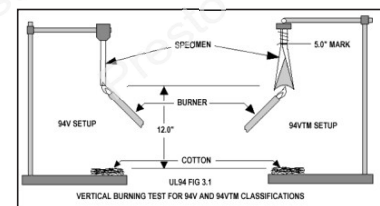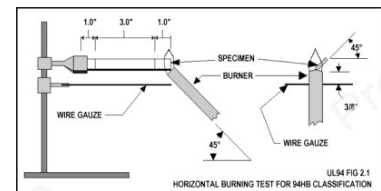


Model No. PFT - UL-94

Flammability Tester UL94 intends this standard to serve as a preliminary indication of a plastics acceptability for use as part of a device or appliance with respect to its flammability. It is not intended to reflect the hazards of a material under actual fire conditions.

UL 94 flammability testing is the first step toward obtaining a plastic recognition and subsequent listing in the "Plastics Recognized Component Directory" (former known as "Yellow Cards"). UL 94 contains the following tests: HB, V-0, V-1, V-2, 5VA or 5VB, VTM-0, VTM-1, or VTM-2, HBF, HBF-1 or HBF-2 .

Burner arrangement is compliant to ASTM D5025. Apparatus has inner volume of 0.5 m3 with exhaust fan.

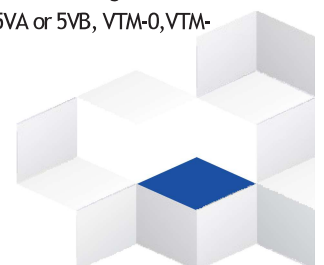



Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India
 P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com



- Faridabad
- Sonipat
- Kolkata
- Mumbai
- Pune
- Ahmedabad
- Chennai
- Bangalore
- Hyderabad

www.prestogroup.com

Key Specifications:

Timer Display	Digital 7 Segment
Distance between Burner & Specimen	Adjustable
Specimen Size	Horizontal Burning Test - 125 x 13 mm Vertical Burning Test - 125 x 13 mm Plaque specimens - 150 x 150 mm Thin material vertical - 200 x 50 mm Foamed horizontal - 150 x 50 mm
Timer accuracy	± 0.5%
Power	220V, Single phase, 50 Hz
Inner dimensions	800x 800 x 800 mm
Digital Timer Range	999.9 Second
Digital Timer Least Count	0.1 second
Igniting Flame Length	Adjustable
Burner angle from horizontal	90 ⁰ , 20 ⁰ , 45 ⁰
Burner nozzle Diameter	9.5 +/- 0.3 mm

Highlights:

Flame Auto Stop	Yes
Material	Stainless Steel
Specimen Holder	1 set
Wire Gauge	Yes
Test type	Horizontal and vertical

- * Horizontal and vertical orientation is possible
- * LPG Empty Cylinder, Regulator, Pipe & Dust Collector Tray not a part of supply

Optional

Flammability Tester with Various Standard DIN 4102-1, IS 15061 also available on request

Thank you customers for choosing us as your partners in growth !



*All the above logos are the sole property of their respective owners and are used purely for depiction purposes only.

OUR OTHER PRODUCTS :-

					
Hot Air Oven Prima	Humidity Chamber	Paint Impact Tester Direct / In-direct	Salt Spray Chamber Korrox - III	Scratch Hardness Tester Motorised	Tensile Testing Machine

