



Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India
P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com

• Faridabad • Sonapat • Kolkata • Mumbai • Pune • Ahmedabad
• Chennai • Bangalore • Hyderabad

www.prestogroup.com



Specific Gravity Balance Digital




Model No. SP 202 (G)

Specific Gravity of any elastomer is of considerable importance for its identification. Besides it can also be used for providing a rough assessment of its physical properties. It can also be used for finding if the sample and the products are prepared from the same compound or not. Displays direct specific Gravity digitally, No calculations required. Based on the Archimedes principle.

Designed for use in laboratory applications, it is easy to use two-button operation, a high-contrast LCD display.

* It is must to sink the sample into the DM water to avoid the floating of sample.

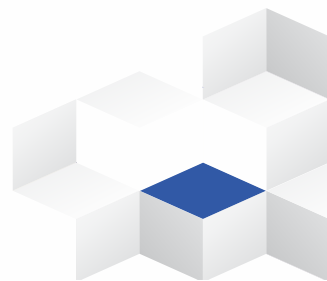



Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India

P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com

• Faridabad • Sonipat • Kolkata • Mumbai • Pune • Ahmedabad
• Chennai • Bangalore • Hyderabad

www.prestogroup.com



Key Specifications:

Max. Weighing Capacity	200 gm
Least Count	0.01 gm
Pan size	120 mm
Display	High contrast LCD display
Calibration	200 gm brass weight (Included)
Linearity	± 0.01 g
Specific Gravity Range	1.0 to 3.0
Medium	Lab Grade DM Water
Power	AC adaptor (Included) or 4 AA battery (not included)

Highlights:

Hanging Hook	Yes
Beaker Platform	Yes
Supply	Main Supply & battery operated

Optional :

- Glass beaker available on request.

Feature :

- Superior shock resistance & overload protection
- Software lock switch & integrated transportation lock
- Adjustable footrest with level indicator.

Thank you customers for choosing us as your partners in growth !



*All the above logos are the sole property of their respective owners and are used purely for depiction purposes only.

OUR OTHER PRODUCTS :-

					
Box Compression Tester	Scuff Resistance Tester	Tensile Tester Digital	Spectrophotometer	Cobb Sizing Tester	Digital GSM Balance

